Made for
 **iPhone**

GENIUS™ 3.0

Receiver-in-Canal Hearing Aid

User Manual

Content

Welcome 4

Your hearing aids 5

Hearing aid type	5
Getting to know your hearing aids	6
Components and names	7
Controls	10
Hearing programs	11
Features	11

Batteries 12

Battery size and handling tips	12
Replacing batteries	13

Daily use 14

Turning on and off	14
Inserting and removing the hearing aids	15
Adjusting the volume	19
Changing the hearing program	19
Further adjustments (optional)	20

Special listening situations 21

On the phone	21
Audio streaming with iPhone	22
Airplane mode	22
Audio induction loops (optional)	23

Maintenance and care 24

Hearing aids	24
Ear pieces	26
Professional maintenance	28

Further information 29

Safety information	29
Accessories	29
Symbols used in this document	29
Troubleshooting	30
Country specific information	32
Service and warranty	33

Welcome

Thank you for choosing Miracle-Ear® hearing aids to accompany you through everyday life. Like anything new, it may take you a little time to become familiar with them.

This guide, along with support from your Hearing Care Specialist, will help you understand the advantages and greater quality of life your hearing aids offer.

To receive the most benefit from your hearing aids, it is recommended that you wear them every day, all day. This will help you get adjusted to them.

CAUTION

It is important to read this user manual and the safety manual thoroughly and completely. Follow the safety information to avoid damage or injury.

Your hearing aids

This user manual describes optional features that your hearing aids may or may not have. Ask your Hearing Care Specialist to indicate the features that are valid for your hearing aids.

Hearing aid type

Your hearing aids are a RIC (Receiver-in-Canal) model. The receiver is placed within the ear canal and connected to the instrument via a receiver cable. The instruments are not intended for children under the age of 3 years or persons with a developmental age of under 3 years.

The wireless functionality enables advanced audiological features and synchronization between your two hearing aids.

Your hearing aids feature *Bluetooth*® low energy* technology allowing for easy data exchange with your smartphone, and for seamless audio streaming with your iPhone.**

* The Bluetooth word mark and logos are owned by the Bluetooth SIG, Inc., and any use of such marks by the legal manufacturer of this product is under licenses. Other trademarks and trade names are those of their respective owners.

** Apple® iPhone is a trademark of Apple Inc., registered in the U.S. and other countries

Getting to know your hearing aids


It is recommended to familiarize yourself with your new hearing aids. With the hearing aids in your hand try using the controls and note their location on the hearing aids. This will make it easier for you to feel and press the controls while wearing the hearing aids.

i If you have problems pressing the controls of your hearing aids while wearing them, you can use a remote control to control your instruments. With a dedicated smartphone app you have even more control options.

Components and names


This user guide describes several types of hearing aids. Use the following pictures to identify the type of hearing aid that you use.

DS RIC 312

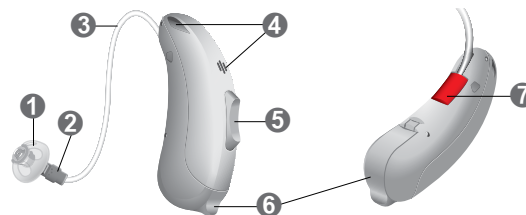- | | |
|---------------------------|---|
| ① Ear piece | ⑥ Battery compartment (on/off switch) |
| ② Receiver | ⑦ Side indicator (red = right ear, blue = left ear) and receiver connection |
| ③ Receiver cable | |
| ④ Microphone openings | |
| ⑤ Rocker switch (control) | |

DS RIC 13

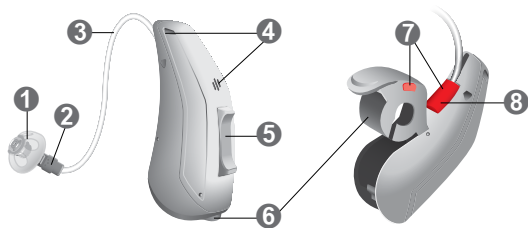- | | |
|----------------------------------|--|
| 1 Ear piece | 6 Battery compartment (on/off switch) |
| 2 Receiver | 7 Side indicator (red = right ear, blue = left ear) |
| 3 Receiver cable | 8 Receiver connection |
| 4 Microphone openings | |
| 5 Rocker switch (control) | |

You can use the following standard ear pieces:


Standard ear pieces		Size
☐		Click Sleeve (vented or closed)
☐		Click Dome™ single (open or closed)
☐		Click Dome semi-open
☐		Click Dome double

You can easily exchange the standard ear pieces. Read more in section "Maintenance and care."

Custom-made ear pieces	
☐	Custom shell
☐	Click Mold™

Controls

With the rocker switch you can adjust the volume or switch hearing programs. Your Hearing Care Specialist has programmed your desired functions to the rocker switch.

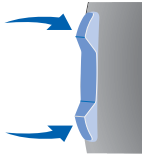Rocker switch function	L	R
Press briefly:		
Program up/down	☐	☐
Volume up/down	☐	☐
Tinnitus control signal level up/down	☐	☐
Press for about 2 seconds:		
Program up/down	☐	☐
Volume up/down	☐	☐
Tinnitus control signal level up/down	☐	☐
Press for more than 3 seconds:		
Standby/turn on	☐	☐
Activate/deactivate airplane mode	☐	☐

L = Left, R = Right

i

You can also use a remote control to change hearing programs and adjust the volume on your hearing aids. With our smartphone app you have even more control options.

Hearing programs

1
2
3
4
5
6

Read more in section “Changing the hearing program.”

Features

- ☐ **Power-on delay** enables whistle-free insertion of hearing aids.
Read more in section “Turning on and off.”
- ☐ The **tinnitus control** feature generates a sound to distract you from your tinnitus.

Batteries

When the battery is low the sound becomes weaker or you will hear an alert signal. The battery type will determine how long you have to replace the battery.

Battery size and handling tips

Ask your Hearing Care Specialist for recommended batteries.

Battery size:

○ 312

○ 13

- Always use the correct battery size for your hearing aids.
- Remove the batteries if you intend not to use the hearing aids for several days.
- Always carry spare batteries.
- Remove empty batteries immediately and follow your local recycling rules for battery disposal.


Replacing batteries

Removing the battery:

- ▶ Open the battery compartment.

- ▶ Use your finger or an optional magnet stick to remove the battery. The magnet stick is available as an accessory from your Hearing Care Specialist.

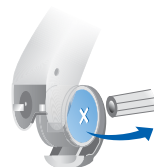Inserting the battery:

- ▶ If the battery has a protective film, remove it only when you are ready to use the battery.
- ▶ Insert the battery with the "+" symbol facing up (refer to the picture).

- ▶ Carefully close the battery compartment. If you feel resistance, the battery is not inserted correctly.

Do not attempt to close the battery compartment by force. It can be damaged.

Daily use

Turning on and off

You have the following options to turn your hearing aids on or off.

Via battery compartment:

- ▶ Turning on: Close the battery compartment.

The **default** volume and hearing program are set.

- ▶ Turning off: Open the battery compartment to the first catch.

Via rocker switch:

- ▶ Turning on or off: Press and hold the rocker switch for several seconds. Refer to section "Controls" for rocker switch settings.

After turning on, the previously used volume and hearing program remain set.

When wearing the hearing aids, an alert tone can indicate when a hearing aid is being turned on or off.


When the **power-on delay** is activated, the hearing aids turn on after a delay of several seconds. During this time you can insert the hearing aids into your ears without experiencing unpleasant feedback whistling.

The "power-on delay" can be activated by your Hearing Care Specialist.

Inserting and removing the hearing aids

Your hearing aids have been fine-tuned for your right and left ear. Colored markers indicate the side:


- Red marker = right ear
- Blue marker = left ear

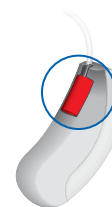
Inserting a hearing aid:

- ▶ For Click Sleeves, make sure that the bend of the Click Sleeve is in line with the bend of the receiver wire.

Correct:

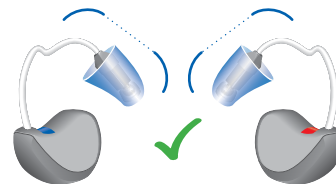Incorrect:

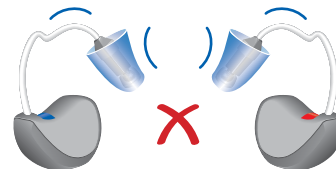
- ▶ Hold the receiver wire at the bend closer to the ear piece.

- ▶ Carefully push the ear piece in the ear canal ①.

- ▶ Twist it slightly until it sits well in the ear canal.

Open and close your mouth to avoid accumulation of air in the ear canal.

- ▶ Lift the hearing aid and place it over the top of your ear ②.

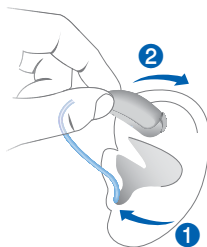CAUTION

Risk of injury!

- ▶ Insert the ear piece carefully and not too deeply into the ear.

- It may be helpful to insert the right hearing aid with the right hand and the left hearing aid with the left hand.
- If you have problems inserting the ear piece, use the other hand to gently pull your earlobe downwards. This opens the ear canal and eases insertion of the ear piece.


The optional retention strand helps to securely retain the ear piece in your ear. To position the retention strand:

- ▶ Bend the retention strand and position it carefully into the lower bowl of your outer ear (refer to the picture).

Removing a hearing aid:

- ▶ Lift the hearing aid and slide it over the top of your ear ①.
- ▶ If your hearing aid is equipped with a custom shell or a Click Mold, remove it by pulling the small removal cord toward the back of your head.
- ▶ For all other ear pieces: Grip the receiver in the ear canal with two fingers and pull it carefully out ②.

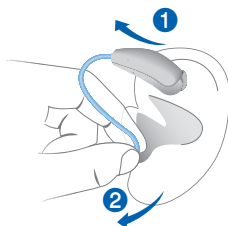Do not pull the receiver cable.

CAUTION

Risk of injury!

- ▶ In very rare cases the ear piece could remain in your ear when removing the hearing aid. If this happens, have the ear piece removed by a medical professional.

Clean and dry your hearing aids after usage. Read more in section "Maintenance and care."

Adjusting the volume

Your hearing aids automatically adjust the volume to your listening situation.

- ▶ If you prefer manual volume adjustment, press the rocker switch.

Refer to section "Controls" for the setting of the rocker switch.

An optional signal tone can indicate the volume change. Once the maximum or minimum volume has been reached, you may hear an optional signal tone.

Changing the hearing program

Depending on the listening situation, your hearing aids automatically adjust their sound.

Your hearing aids may also have several hearing programs which allow you to change the sound, if needed. An optional signal tone can indicate the program change.

- ▶ To change the hearing program, press the rocker switch.

Refer to section "Controls" for the setting of the rocker switch. Refer to section "Hearing programs" for a list of your hearing programs.

Further adjustments (optional)


You can also use the hearing aid controls to change other features, for example, the level of the tinnitus therapy signal.

Refer to section "Controls" for the setting of your controls.

Special listening situations

On the phone

When you are on the phone, hold the telephone receiver slightly above your ear. The hearing aid and telephone receiver have to be aligned. Turn the receiver slightly so that the ear is not completely covered.

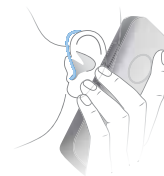Telephone program

You may prefer a certain volume when using the phone. Ask your Hearing Care Specialist to configure a telephone program.

- Switch to the telephone program whenever you are on the phone.

If a telephone program is configured for your hearing aids, it is listed in section "Hearing programs."

Audio streaming with iPhone

Your hearing aids are **Made for iPhone®** hearing aids. This means you can get phone calls and listen to music from your iPhone directly to your hearing aids.

For more information about compatible iOS devices, pairing, audio streaming and other useful functions, contact your Hearing Care Specialist.

Airplane mode

In areas where the use of Bluetooth technology is not allowed (e.g. on some airplanes), you can activate airplane mode. This turns off the *Bluetooth* in your hearing aids temporarily. The hearing aids still work without Bluetooth, however direct audio streaming will not be possible and other functions will not be available.

- To activate or deactivate the airplane mode, use the MEcontrol smartphone app.

An alert tone indicates when airplane mode is activated or deactivated.

Audio induction loops (optional)

Some telephones and public places like theatres offer the audio signal (music and speech) via an audio induction loop. With this system, your hearing aids can directly receive the desired signal – without distracting noise in the environment.

You can usually recognize audio induction loop systems by certain signs.

If you want to use this system, a special battery compartment is required. It is available as an accessory only for certain types of hearing aids. Consult your Hearing Care Specialist for more information.

Whenever you are in a place with an audio induction loop:

- Switch to the telecoil program.

If a telecoil program is configured for your hearing aids, it is listed in section "Hearing programs".

Telecoil program and remote control

NOTICE

For hearing aids with wireless connectivity: When the telecoil program is active, the remote control may cause a pulsing noise.

- Use the remote control with a distance greater than 10 cm (4 inches).

Maintenance and care

To prevent damage it is important that you take care of your hearing aids and follow a few basic instructions, which will soon become a part of your daily routine.

Hearing Aids

Drying and storage

- ▶ Dry your hearing aids overnight.
- ▶ Ask your Hearing Care Specialist for recommended drying products.
- ▶ For longer periods of non-use, store your hearing aids with an open battery compartment and place in a drying system to avoid the adverse effects of moisture.

- ▶ Never use running water or immerse the devices in water.
- ▶ Never apply pressure while cleaning.
- ▶ Ask your Hearing Care Specialist for recommended cleaning products, special care sets, or more information on how to keep your hearing aids in good condition.

Ear pieces


Cleaning

Cerumen (ear wax) may accumulate on the ear pieces and may affect sound quality. Clean the ear pieces daily.

For all types of ear pieces:


- Clean the ear piece right after removal with a soft and dry tissue.

This prevents cerumen from becoming dry and hard, causing build-up on the ear piece.

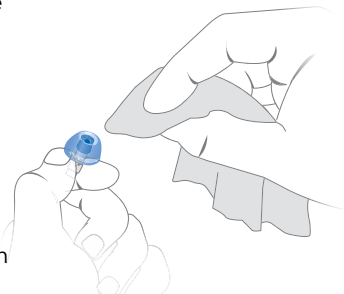For Click Domes or Click Sleeves:

- Squeeze the tip of the ear piece.

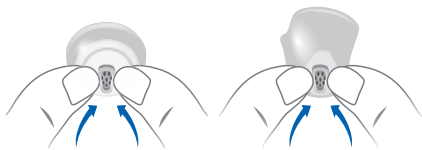
Exchanging

Exchange standard ear pieces approximately every three months. Replace them sooner if you notice cracks or other changes. The procedure for exchanging standard ear pieces depends on the type of ear piece. In section "Components and names", your Hearing Care Specialist has marked your type of ear piece.

Exchanging Click Domes

- To remove the old Click Dome, use the removal tool or follow the instructions on the packaging of the Click Domes. Be careful not to pull the receiver cable.

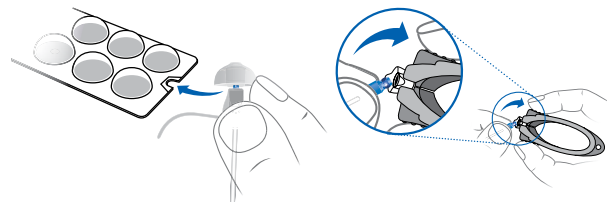The removal tool is available as an accessory.

- Make sure that the new Click Dome clicks noticeably in place.
- Check to make sure that the click dome has a firm fit.

Exchanging Click Sleeves

- ▶ Turn the Click Sleeve inside out.
- ▶ Grasp the Click Sleeve and the receiver (not the receiver wire) with your fingers and pull the Click Sleeve off.
- ▶ Make sure that the new Click Sleeve clicks noticeably in place.
- ▶ Check to make sure that the click sleeve has a firm fit.

Professional maintenance

Your Hearing Care Specialist can perform a thorough cleaning and maintenance.

Custom-made ear pieces and wax filters should be exchanged as needed by your Hearing Care Specialist.

Ask your Hearing Care Specialist for a recommendation of maintenance intervals and support.

Further information

Safety information

For more safety information, refer to the safety manual supplied with the hearing aids.

Accessories

You can use the MEcontrol App for convenient remote control of your hearing aids. The app also allows you to deactivate signal tones, to create individual sound profiles and many more functions.

You can also use a remote control.

Consult your Hearing Care Specialist for further information.

Symbols used in this document

Points out a situation that could lead to serious, moderate, or minor injuries.

Indicates possible property damage.

Advice and tips on how to handle your device better.

"Made for iPhone" means that an electronic accessory has been designed to connect specifically to an iPhone and has been certified by the developer

to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Please note that the use of this accessory with iPhone may affect wireless performance.

Troubleshooting

Problem and possible solutions

Sound is weak.

- Increase the volume.
 - Replace the empty battery.
 - Clean or exchange the ear piece.
-

Hearing aid emits whistling sounds.

- Reinsert the ear piece until it fits securely.
 - Decrease the volume.
 - Clean or exchange the ear piece.
-

Sound is distorted.

- Decrease the volume.
 - Replace the empty battery.
 - Clean or exchange the ear piece.
-

Hearing aid emits signal tones.

- Replace the empty battery.
-

Problem and possible solutions

Hearing aid does not work.

- Turn the hearing aid on.
 - Gently close the battery compartment completely.
 - Replace the empty battery.
 - Make sure the battery is inserted correctly.
 - Power-on delay is active. Wait several seconds and check again.
-

Consult your Hearing Care Specialist if you encounter further problems.

Country Specific Information

USA and Canada

Certification information for hearing aids with wireless functionality:

GENIUS™ 3.0 DS RIC 312

- FCC ID: SGI-RIC002
- IC: 267AB-RIC002

GENIUS™ 3.0 DS RIC 13 (without telecoil option)

- FCC ID: SGI-RIC001
- IC: 267AB-RIC001

GENIUS™ 3.0 DS RIC 13 (with telecoil option)

- Contains FCC ID: SGI-RFM001, SGI-MIM001
- IC: 267AB-WP4N1

Service and warranty

Serial numbers	Hearing aid	Receiver unit
Left:		
Right:		

Service dates	
1:	4:
2:	5:
3:	6:

Warranty
Date of purchase:
Warranty period in months:

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Copyright ©2019 Miracle-Ear, Inc.
Distributed by Miracle-Ear, Inc.
Derechos de autor ©2019 Miracle-Ear, Inc.