MEENERGY™ IN Sync

User Manual

Content

Welcome	4
Your hearing aids	5
Hearing aid type	5
Getting to know your hearing aids	6
Components and names	6
Controls	6
Hearing programs	8
Features	8
Daily use	9
Charging	9
Turning on and off hearing aids	10
Switching to standby mode	11
Inserting and removing the hearing aids	12
Adjusting the volume	14
Changing the hearing program	15
Further adjustments (optional)	15
Special listening situations	16
On the phone	16
Audio streaming with iPhone	17
Audio streaming with Android phone	17
Airplane mode for hearing aids	17

Maintenance and care	18
Cleaning	18
Drying	18
Storage	19
Professional maintenance	20
Important safety information	21
Safety information on power cells	21
Further information	24
Accessories	24
Operating, transport and storage conditions	24
Explanation of symbols	26
Troubleshooting	28
Country-specific information	30
Service and warranty	31

Welcome

Thank you for choosing Miracle-Ear® hearing aids to accompany you through everyday life. Like anything new, it may take you a little time to become familiar with them.

This guide, along with support from your Hearing Care Professional, will help you understand the advantages and greater quality of life your hearing aids offer.

To receive the most benefit from your hearing aids, it is recommended that you wear them every day, all day. This will help you get adjusted to them.

The device may not look exactly like the guiding illustrations in these instructions. We reserve the right to make any changes we consider necessary.

CAUTION

It is important to read this user manual and the safety manual thoroughly and completely. Follow the safety information to avoid damage or injury.

Before using your hearing aids for the first time, fully charge them.

Your hearing aids

This user manual describes optional features that your hearing aids may or may not have.

Ask your Hearing Care Professional to indicate the features that are valid for your hearing aids.

Hearing aid type

Your hearing aids are worn in the outer ear. The hearing aids are not intended for children under the age of 3 years or persons with a developmental age of under 3 years.

A power cell (lithium-ion rechargeable battery) is permanently built into your hearing aid. This enables easy charging with your charger.

The wireless functionality enables advanced audiological features and synchronization between your two hearing aids.

Your hearing aids feature Bluetooth® Low Energy* technology allowing for easy data exchange with your smartphone, and for seamless audio streaming with your iPhone.**

* The Bluetooth word mark and logos are owned by the Bluetooth SIG, Inc., and any use of such marks by WS Audiology Denmark A/S is under license. Other trademarks and trade names are those of their respective owners.


** iPad, iPhone, and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries.

Getting to know your hearing aids

Please download the the Miracle-Ear® App for remote control features.

You can also ask your Hearing Care Professional about another remote control option.

Components and names

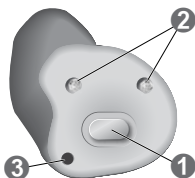❶ Push button

❸ Vent

❷ Microphone opening

Controls

With the push button you can, for example, switch hearing programs. The Hearing Care Professional has programmed your desired functions to the push button.

You can use the Miracle-Ear® App for remote control.

Push button function	L	R
Press briefly:		
Program up	○	○
Volume up	○	○
Volume down	○	○
Tinnitus therapy signal level up	○	○
Tinnitus therapy signal level down	○	○
TV stream on/off	○	○
Press for about 2 seconds:		
Program up	○	○
Volume up	○	○
Volume down	○	○
Tinnitus therapy signal level up	○	○
Tinnitus therapy signal level down	○	○
TV stream on/off	○	○
Activate/deactivate airplane mode	○	○
Press for more than 3 seconds:		
Standby/turn on	○	○
L = Left, R = Right		

Hearing programs

1

2

3

4

5

6

Read more in section "Changing the hearing program".

Features


- ☐ The tinnitus therapy feature generates a sound to distract you from your tinnitus.

Daily use

Charging

Use MECHARGE™ compact to charge your hearing aids. Place your hearing aids in the charger.

Follow the instructions in the charger's user manual. You can also find useful tips for charging in the charger's user manual.

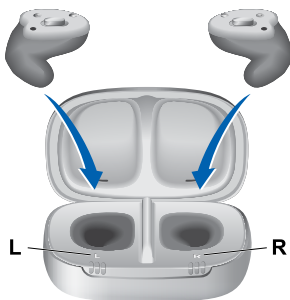Turning on and off hearing aids

You have the following options to turn your hearing aids on or off.

► Turning on:

Take the hearing aids out of the charger.

The startup melody is played in your hearing aids.

The default volume and hearing program are set.

► Turning off in charger:

Place the hearing aids in the charger.


Make sure the charger is on

- If the charger is turned off and you insert hearing aids that are turned on, these are not automatically turned off and will not be recharged.
- If the charger is on and you insert hearing aids that are turned on, these are automatically turned off and recharged.

► Turning off on push button:

Use the push button on the hearing aid.

Simply press and hold for 10 seconds.

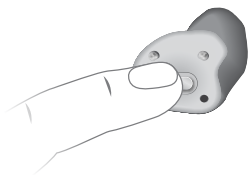Switching to standby mode

Via remote control, you can switch your hearing aids to standby mode. This mutes the hearing aids. When you leave standby mode, the previously used volume and hearing program are set.

- ▶ In standby mode, the hearing aids are not turned off completely. They draw some power. Therefore we recommend to use standby mode only for a short period of time.
- ▶ If you want to leave the standby mode but the remote control is not at hand:

Turn your hearing aids off and on again by shortly placing them in the charger until one or more LEDs light up.

To do this, the charger must be powered. Wait until the hearing aids play the startup melody. This may take several seconds. Note that the default volume and hearing program are set.

Inserting and removing the hearing aids


Your hearing aids have been fine-tuned for your right and left ear. The labeling on the hearing aids indicates the side:

- "R" = right ear
- "L" = left ear

Inserting a hearing aid:

- ▶ Carefully push the hearing aid into the ear canal.
- ▶ Twist it slightly until it sits well.


Open and close your mouth to avoid accumulation of air in the ear canal.

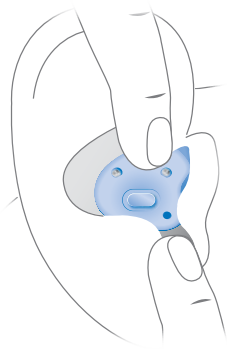Removing a hearing aid:

- ▶ Push lightly on the back of your ear to loosen the hearing aid.

- Grip the hearing aid in the ear canal with two fingers and pull it carefully out.

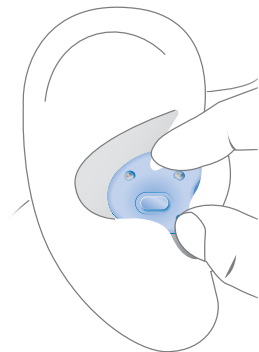Clean and dry your hearing aids after usage. Read more in section "Maintenance and care".

Adjusting the volume

Your hearing aids automatically adjust the volume to the listening situation.

- ▶ If you prefer manual volume adjustment, use a remote control option.

An optional signal can indicate the volume change.

Changing the hearing program

Depending on the listening situation, your hearing aids automatically adjust their sound.

Your hearing aids may also have several hearing programs which allow you to change the sound, if needed. An optional signal tone can indicate the program change.

- To change the hearing program manually, use a remote control option.

Refer to section "Hearing programs" for a list of your hearing programs.

The default volume is applied.

Further adjustments (optional)


You can also use the hearing aid controls to change other features, for example, the level of the tinnitus therapy signal.

Refer to section "Controls" for the setting of your controls.

Special listening situations

On the phone

When you are on the phone, turn the receiver slightly so that it does not completely cover your ear.

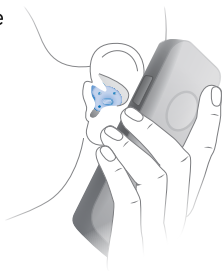Telephone program

You may prefer a certain volume when using the phone. Ask your Hearing Care Professional to configure a telephone program.

- Switch to the telephone program whenever you are on the phone.

If a telephone program is configured for your hearing aids, it is listed in section "Hearing programs".

Audio streaming with iPhone

Your hearing aids are Made for iPhone hearing aids. This means you can get phone calls and listen to music from your iPhone directly into your hearing aids.

For more information about compatible iOS devices, pairing, audio streaming and other useful functions, contact your Hearing Care Professional.

Audio streaming with Android phone

If your mobile phone supports Audio Streaming for Hearing Aids (ASHA), you can get phone calls and listen to music from your mobile phone directly into your hearing aids.

For more information about compatible Android devices, pairing, audio streaming and other useful functions, contact your Hearing Care Professional.

Airplane mode for hearing aids

Airplane mode is when Bluetooth® of your hearing aids is turned off. In areas where the use of Bluetooth technology is not allowed (e.g. on some airplanes), you can deactivate Bluetooth mode. This turns Bluetooth in your hearing aids off temporarily. The hearing aids still work without Bluetooth, however direct audio streaming will not be possible and other functions will not be available.

- To deactivate or activate the Bluetooth mode, use the smartphone app.

An alert tone indicates when Bluetooth mode is deactivated or activated.


Maintenance and care

To prevent damage it is important that you take care of your hearing aids and follow a few basic rules, which will soon become a part of your daily routine.

Cleaning

For hygiene reasons and to maintain functionality, clean your hearing aids daily.

- ▶ Clean your hearing aids daily with a soft, dry tissue.
- ▶ Your hearing aids may have a wax protection. This prevents the ingress of earwax. Be careful not to damage the wax protection when cleaning your hearing aids.

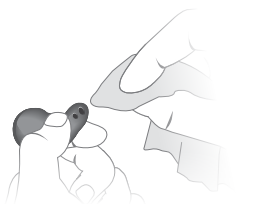- ▶ Never use running water or immerse the devices in water.
- ▶ Never apply pressure while cleaning.
- ▶ Ask your Hearing Care Professional for recommended cleaning products, special care sets, or more information on how to keep your hearing aids in good condition.

Drying

- ▶ Dry your hearing aids overnight.
- ▶ Ask your Hearing Care Professional for more information.

Storage

Short-term storage (up to several days):

- ▶ Turn off your hearing aids by placing them in the charger.

The charger must be on. If the charger is not on and you place your hearing aids in the charger, the hearing aids are not turned off.

- ▶ Note that when you turn off the hearing aids via a remote control or smartphone app, the hearing aids are not turned off completely. They are in standby-mode and still draw some power.

Long-term storage (weeks, months, ...):

- ▶ Fully charge your hearing aids; then leave them in the charger with the lid closed. Once the hearing aids are fully charged, the charger and the hearing aids are turned off automatically.
- ▶ If you need to store your hearing aids for more than 3 months, put them in the charger and make sure the charger is turned off.
- ▶ We recommend using drying products while storing the hearing aids.
- ▶ For longer periods of non-use, store your hearing aids in a dry environment.
- ▶ Observe the storage conditions in section "Operating, transport and storage conditions".

Professional maintenance

Your Hearing Care Professional can perform a thorough professional cleaning and maintenance.

In the event that the wax protection or microphone protection (windscreen) needs to be replaced, this should be carried out by a professional.

Ask you Hearing Care Professional for your individual recommendation of maintenance intervals and for support.

Important safety information

WARNING

Risk of strangulation!

- ▶ Ensure that the lanyard does not get entangled, e. g. in moving mechanical parts of machines.
- ▶ Use only the delivered breakaway lanyard. It will automatically detach when entangled.
- ▶ Do not use the lanyard for children under the age of twelve.

Safety information on power cells

This section contains important safety information regarding the power cell. For more safety information, refer to the safety manual supplied with the hearing aids.

WARNING

A power cell (lithium-ion rechargeable battery) is permanently built into your hearing aid. Incorrect use of the power cell or the hearing aid can cause the power cell to burst.

Risk of injury, fire or explosion!

- ▶ Follow the safety instructions for power cells in this section.
- ▶ Observe the operating conditions.
- ▶ Pay special attention to protect your hearing aids from extreme heat (e.g. fire, microwave oven, inductive oven or other high-induction fields) and sunlight.
- ▶ Charge your hearing aids only with the approved charger. Ask your Hearing Care Professional for support.
- ▶ Do not disassemble the power cell.

WARNING

Incorrect use of the power cell or the hearing aid can cause the power cell to burst.

Risk of injury, fire or explosion!

Replacement and repair have to be performed by Service Center only.

- ▶ Do not repair or replace the power cell yourself.
- ▶ Do not use the hearing aids when their performance after charging reduces considerably.
- ▶ In the unlikely event of a power cell bursting or exploding, keep all combustible materials away.

NOTICE

- ▶ To avoid environmental pollution, do not throw power cells, hearing aids or charger into household trash.
- ▶ Recycle or dispose of power cells, hearing aids or charger according to local regulations.

For more safety information than included in this user manual, refer to the safety manual supplied with the device.

Further information

Accessories

The type of available accessories depends on the hearing aid type. Examples are: a remote control, an audio streaming device, or a smartphone app for remote control.

Consult your Hearing Care Professional for further information.

Operating, transport and storage conditions

Note that you have to charge the hearing aids at least every 6 months. For details refer to section "Maintenance and care".

Operating conditions	Discharging	Charging
Temperature	0-50 °C 32-122 °F	10-35 °C 50-95 °F
Relative humidity	5-93 %	5-93 %
Atmospheric pressure	700-1060 hPa	700-1060 hPa

During extended periods of transport and storage, please observe the conditions below.

	Storage	Transport
Temperature (recommended)	10-25 °C 50-77 °F	-20-60 °C -4-140 °F
Temperature (maximum range)	10-40 °C 50-104 °F	-20-60 °C -4-140 °F
Relative humidity (recommended)	20-80 %	5-90 %
Relative humidity (maximum range)	10-80 %	5-90 %
Atmospheric pressure	700-1060 hPa	700-1060 hPa

WARNING

Risk of injury, fire or explosion during shipment.

- Ship power cells, hearing aids and/or charger according to local regulations.

Explanation of symbols

Symbols used in this document

 Points out a situation that could lead to serious, moderate, or minor injuries.

 Indicates possible property damage.

 Advice and tips on how to handle your device better.

 Frequently asked question.

Symbols on the device or packaging

 Indicates the legal manufacturer of the device.

 Do not dispose of the device with general domestic waste. Read more in section "Disposal information".

 Packaging can be recycled.

 Read and follow the instructions in the user manual.

 Keep dry during transportation.

 Temperature during transportation.
Read more in section "Operating, transport and storage conditions".

Symbols on the device or packaging

Relative humidity during transportation. Read more in section "Operating, transport and storage conditions".

Atmospheric pressure during transportation. Read more in section "Operating, transport and storage conditions".

Medical device.

Bluetooth® Low Energy technology.

Made for

iPhone | iPad | iPod

"Made for iPhone", "Made for iPad", and "Made for iPod" mean that an electronic accessory has been designed to connect specifically to iPhone, iPad, or iPod, respectively, and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Please note that the use of this accessory with iPhone, iPad, or iPod may affect wireless performance.

Charger data

Under the removable silicone lining in the lid of the charger is a label with symbols, serial number, CE mark and legal manufacturer of the charger. Lift the lining from the lid with your fingernail and pull it off. To reinsert, press around the edges.

Troubleshooting

Problem and possible solutions

Sound is weak.

- Increase the volume.
 - Charge the hearing aid.
 - Clean or exchange the ear piece.
-

Hearing aid emits whistling sounds.

- Reinsert the ear piece until it fits securely.
 - Decrease the volume.
 - Clean or exchange the ear piece.
-

Sound is distorted.

- Decrease the volume.
 - Charge the hearing aid.
 - Clean or exchange the ear piece.
-

Hearing aid emits signal tones.

- Charge the hearing aid.
-

Problem and possible solutions

Hearing aid does not work.

- Turn the hearing aid on.
 - Charge the hearing aid. Observe also the status LEDs on the charger.
-

Consult your Hearing Care Professional if you encounter further problems.

Country-specific information

USA and Canada

Certification information for hearing aids with wireless functionality:

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Copyright ©2022 Miracle-Ear, Inc.
Distributed By Miracle-Ear, Inc.
Derechos de autor ©2022 Miracle-Ear, Inc.